

Illustration Instruction
MD680st/MD660 series
Sealers

MD680st,MD680st WA,MD680st WM,MD660
封口机图示说明

J&D Group Germany
www.md-ca.com

Medicare Co. Ltd. Shenzhen
深圳市美康达医疗保健有限公司
www.mkd1995.com

MD680st/MD660 series

MD680st WA and MD680st WM refer to them.



MD680st


MD660



MD680st/MD660 series MD680st WA and MD680st WM refer to them.

Replace the heater tapes.
更换加热条胶带。


Dismantle 4 screws and open the cover. 将前后四颗螺丝取下并打开上盖。

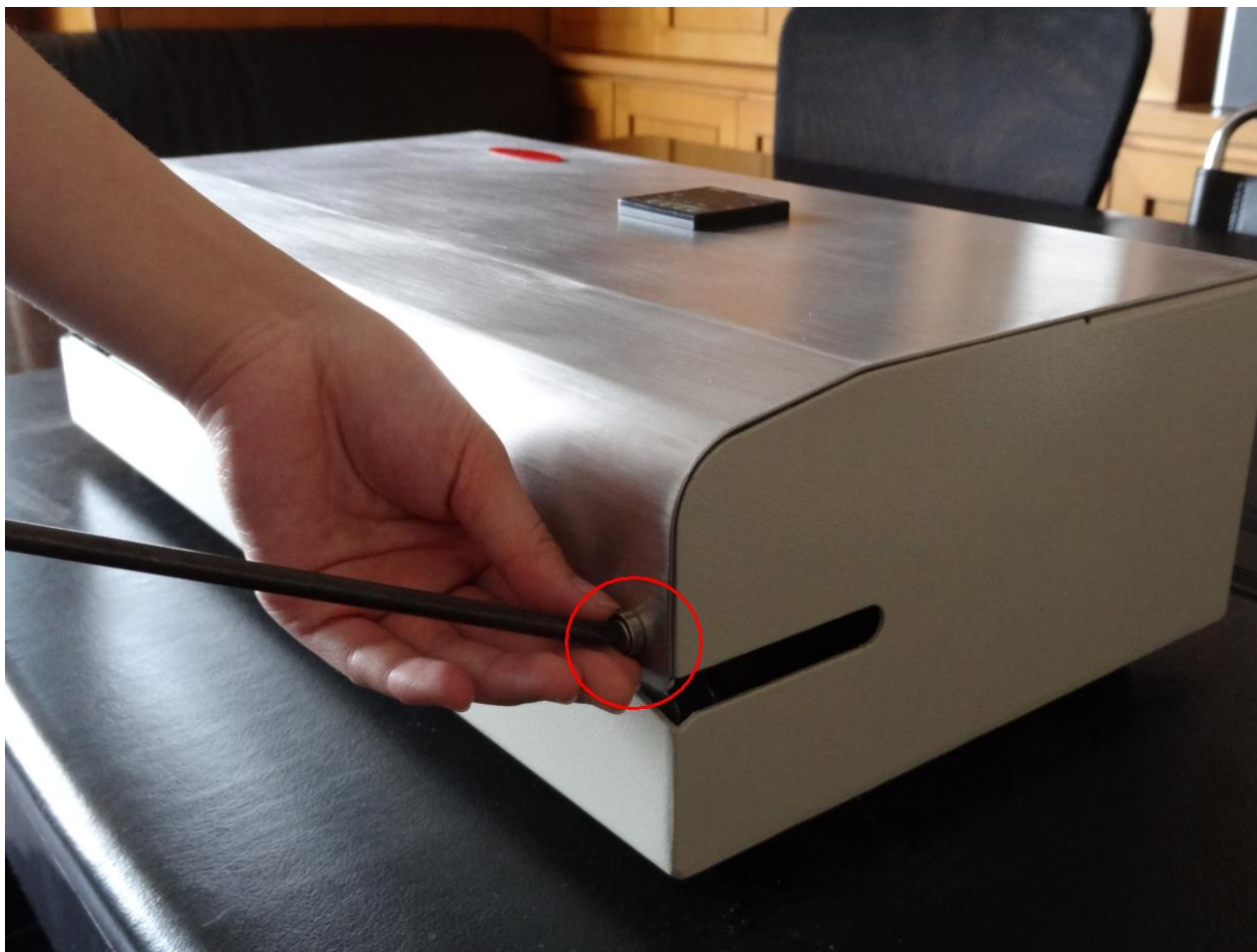

Dismantle 2 screws at the base house.
卸下底座两个螺丝。



Dismantle 3 screws and earth line, see the red circles. 将3枚固定螺丝与地线（见红圈）卸下。

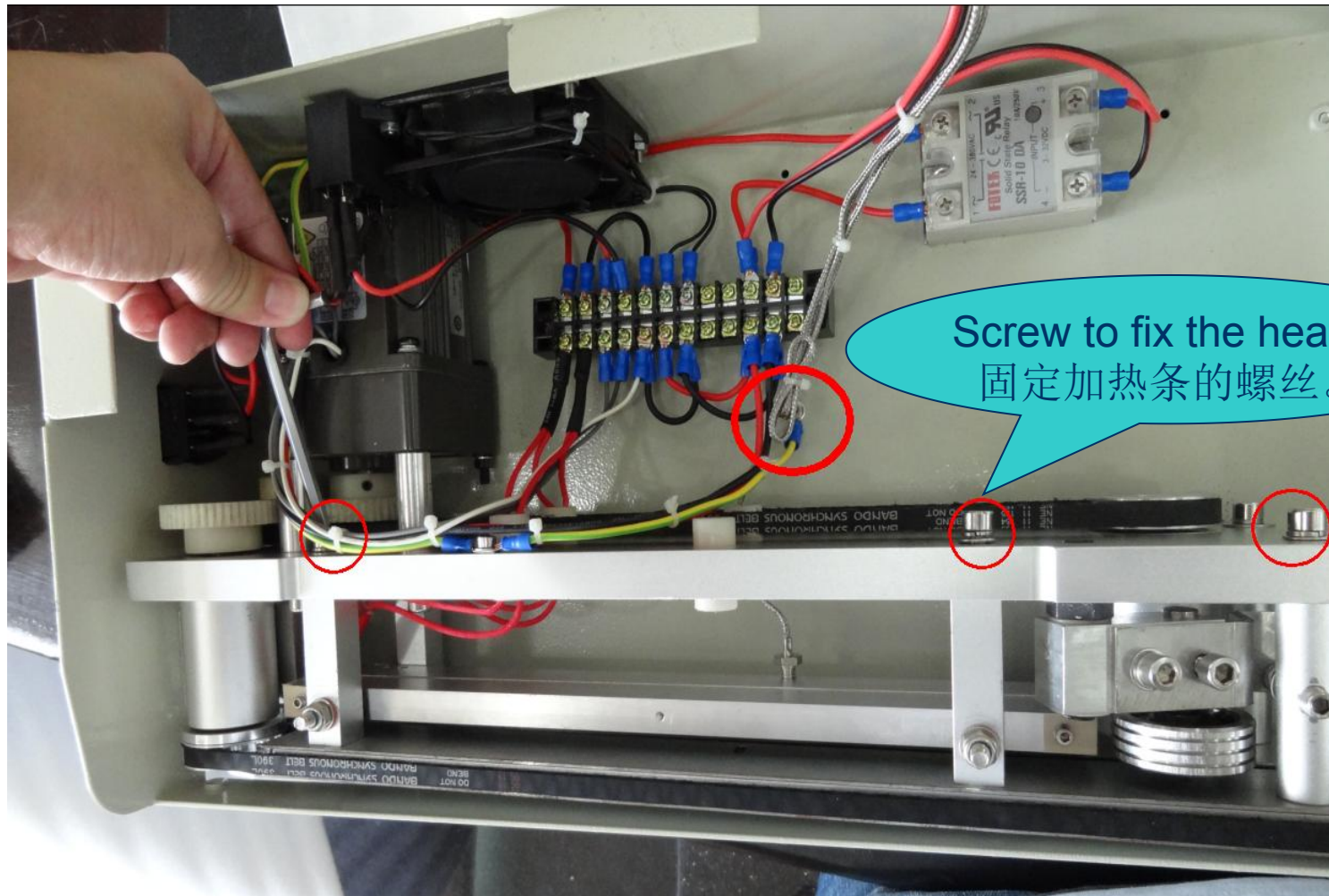

Dismantle the heaters. 将加热条卸下。



Replace the tapes. 更换胶带。

- Take out the tape from the heaters, clean the heater and then replace a new tape.
 - The same way to replace the tapes of upper and lower heaters.
 - Reverse the above steps and install the sealer.
-
- 从加热条上撕下胶带，清洁加热条，然后换上新的胶带。
 - 更换上下加热条的胶带方法相同。
 - 按上述倒回步骤安装好封口机。



**Call your distributor for more
information**

Thanks

如有需要，请联系您的经销商，谢谢！

V3.0 2024